

Linda3 chair embodies the symbol of Tecnodent's philosophy. Enduring elegance and design, soft pure lines, perfect combination of polyurethane and metal for absolute reliability and optimal ergonomics.

Linda3 provides total liberty of movements, also for little rooms. Anatomical and ultra-slim backrest, with lateral armrest support wings, ensures total comfort to the patient and easy access to the chair for the professional. Folding leg rest and the minimum height of 51 cm guarantee easiest access also to elder or disabled patients.

Seamless high quality leather-like fabric provides reassuring comfort and quick and reliable hygiene.

Joysticks on the base allow direct and easy control of chair functions and a quick recall of memorized working positions.

Linda3 can also be equipped with rotating base, in case of further diagnostic devices use.

STANDARD AND OPTIONAL CHARACTERISTICS

• Standard • Optional

	2009new	Lux it	ECO.next	Linda3	Sting
Headrests	ELLE FISSO	•		•	
	ELLE/ELLE2 MOBILE	•	•	•	•
	UNI/Light	•	•	•	•
Pillows	C95	•	•	•	•
	C95 Gel Appeal	•	•	•	•
	C2002	•	•	•	•
Backrest	Narrow	•	•		
	Large	•	•	•	•
	"T" shaped	•			
Armrest	Right armrest	•	•	•	•
	Left armrest	•	•	•	•
Programs	Return to "zero" program	•	•	•	•
	Rinse program	•	•	•	•
Unit connection	Double side drilled upper body	•	•	•	•
	Fixed central arm		•	•	
	Rotation	•	•	•	•
	Rotation with central arm	•	•	•	•
Function	Trendelenburg	•	•	•	•
	Knee-break legrest				•
	Independent legrest movement (3rd motor)				•
	Securities	•	•	•	•
Control	Small integrated utility box		•		
	Large integrated utility box	•		•	•
	Joystick		•	•	•
Control	Foot control on the base	•			
	Mobile foot control	•	•	•	•
	Infrared remote control		•	•	•

HEADRESTS

Tecnodent manual Headrests are available in 4 different models: ELLE FISSO, ELLE/ELLE2 MOBILE, UNI and Light. Every headrest was made for meeting every operator's need.

ELLE/ELLE2. ELLE/ELLE2 FISSO has only the movement along its vertical axis. ELLE/ELLE2 MOBILE is an anatomical headrest, characterized by movement of inclination of the upper part (double articulation) for a correct cervical support.



ELLE/ELLE2



C95



C2002

Magnetic pillows. ELLE/ELLE2 type, both fix and mobile version, can be combined with magnetic pillow C95 or C2002. Either pillow can be applied directly on the backrest, for ensure a better comfort, independently of the patient's height.



UNI



Light


UNI and Light. The UNI and Light headrests are multi-articulated and allows positioning of the headrest along 3 axis.

Magnetic Gel Appeal Pillow adapts to the physical structure of every person, absorbs and distributes head and instrument pressure over the entire contact surface. Long and pleasant sensation of freshness during treatment.



C95 Gel Appeal



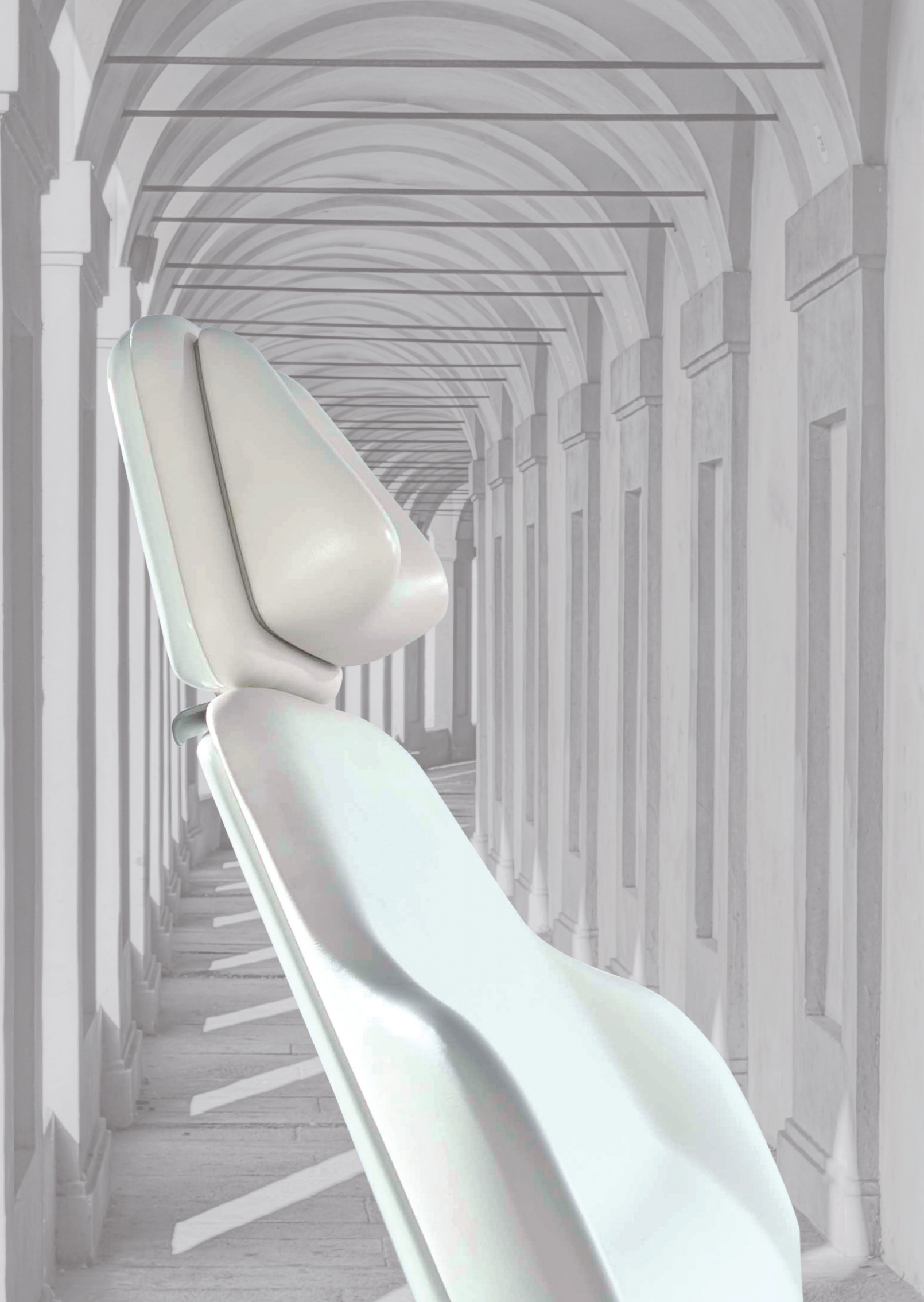

LINDA3

- Electromechanical lifting, lowering and backrest tilting movements synchronized with the leg rest;
- Compensated backrest tilting, up to the horizontal position;
- Anatomic ultra-slim backrest with lateral arm-support wings, for perfect comfort for the patient and liberty of movement for the practitioner;
- Anatomic headrest ELLE FISSO ensures the maximum comfort for the patient and the perfect visibility for the professional;
- Wrapping seamless upholstery in high quality artificial leather;
- Possible customization with different colors of polyurethane and fabric;
- Complete symmetry for right-handed or left-handed professionals;
- Double Joystick control on the base for quick and efficient chair movements with hands free, allowing the assistant perfect control of the chair;

- Security STOPS on the backrest, on the base and on the pantograph arm;
- Return to “zero” position automatic program;
- Rinse position automatic program;
- Further three memorizable positions (OPTION);
- Possibility to equip the chair with one/two rotating armrests to improve the accessibility to the chair and to ensure a natural support for the arms, allowing a relaxed position during the treatment (OPTION);
- RAL 9002 standard chair painting;
- Special RAL painting (OPTION);
- Memory Foam or Hand Made Soft Foam upholstery (OPTION).



TECHNICAL FEATURES

Minimum height	mm	510
Maximum height	mm	840
Backrest tilting angle	°	-5/+70
Rotation angle	°	325
Width	mm	660
Length (without headrest)	mm	1520
Length (with headrest)	mm	1900
Gross weight	Kg	145
Net weight	Kg	115
Maximum lifting capability (patient + unit)	Kg	135+75
Insulation classification		I
Applied parts class		B

LINDA3 - 3Motors

- Independent legrest movement with third motor;
- Leg section at 78° for easy patient access;
- Short Legrest;
- 3 programs + zero and rinse position;
- For all other characteristics and options please refer to standard Linda3



Right/left armrest (OPTION)

Rotating armrests to ease patient's access to the chair.



Standard base

Compact polyurethane base which guarantees stability and liberty of access and movements to the professional.



“C” Metal base (OPTION)


Specific compact metal base (RAL7016) which enable easy access to the chair.
Possibility of different painting (Extra Option)

Swivel base (OPTION)

Allows 325° rotation around the vertical axis.

Infrared control (OPTION)

Manual control of movements and memorized positions of the chair.
(not possible for Linda3 3 motors)

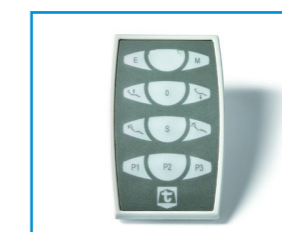

Mobile foot control “Easy Touch” (OPTION)

Easy and intuitive foot control which allows the professional precise hands-free positioning of the chair, for the maximum hygiene.